

ENDOCOR
PIONEERING INNOVATION FOR LIFE™

ENDOGLIDE™

PTA BALLOON CATHETER



LOWEST TIP ENTRY
PROFILE IN THE WORLD

0.035" 0.018" 0.014"

WITH AN IMPRESSIVE RANGE OF SIZES

HIGHLIGHTS

LOWEST TIP ENTRY PROFILE:
THROUGH OUR SPECIAL BALLOON TECHNOLOGY ENDOGLIDE™ ACHIEVES ULTRA LOW TIP ENTRY PROFILE AND ENSURES CROSSING THE MOST DIFFICULT LESIONS EASILY.

WIDEST RANGE OF SIZES:
ALMOST LIMITLESS COMBINATION OPTIONS WITH 0.035" / 0.018" / 0.014" GUIDEWIRE COMPATIBILITY AND AVAILABLE BALLOON LENGTHS FROM 10MM TO 280MM AND DIAMETERS FROM 1.25MM TO 12MM GUARANTEE BEST PATIENT TREATMENT.

MAXIMUM FLEXIBILITY:
ENDOGLIDE™ EXCELS WITH ITS ULTRA FLEXIBILITY AND ENSURES EASY NAVIGATION THROUGH VERY TORTUOUS ANATOMY TO REACH THE MOST DISTAL LESIONS.

FAST INFLATION AND DEFLATION TIMES:
MAXIMISES EFFICIENCY DURING TREATMENT AND REDUCES THE PROCEDURE TIME COMPARED TO OTHER INDUSTRY PARTICIPANTS.

LESS RADIATION EXPOSURE:
REDUCES RADIATION LOAD FOR CLINICIAN, TEAM AND PATIENT DUE TO DECREASED TIME FOR BALLOON OPERATION.

ADVANCED PUSHABILITY:
PUSHABILITY IS INCREASED FOR A SIGNIFICANT AMOUNT OF APPLIED FORCE TRANSITION FROM THE HUB TO THE TIP.

HIGH RATED BURST PRESSURE:
OPTIMAL ANGIOGRAPHICAL RESULTS AND INCREASED SAFETY LEVEL FOR THE PATIENTS IS ENSURED THROUGH HIGH RATED BURST PRESSURE.




TIP COMPARISON OF 0.014" COMPATIBLE PRODUCTS

| | ABBOTT | MEDTRONIC | CORDIS | EV3 | ENDOCOR |
|---------------------|--------|-----------|--------|-----|---------|
| TIP | | | | | |
| TIP ENTRANCE ø (µm) | 606 | 456 | 642 | 404 | 395 |

ENDOGLIDE™ OFFERS THE LOWEST TIP ENTRANCE DIAMETER OF 0.014" COMPATIBLE PRODUCTS IN THE WORLD

TIP ENTRANCE DIAMETER OF 0.018" COMPATIBLE PRODUCTS

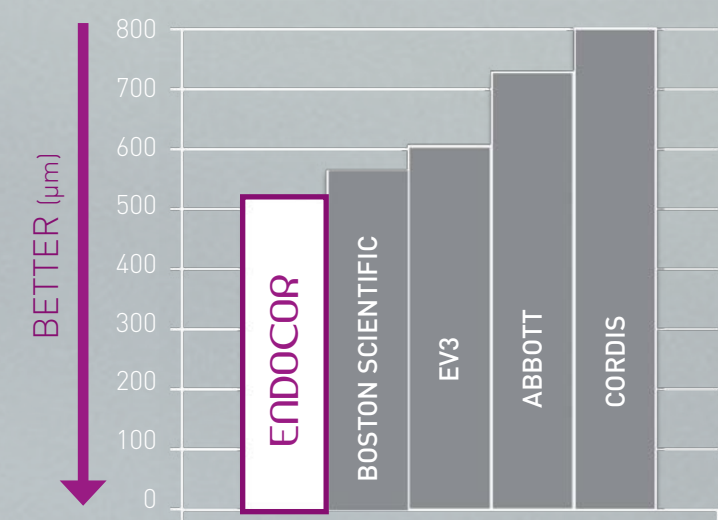

ENDOGLIDE™ ALSO EXCELS IN 0.018" COMPATIBLE PRODUCTS AGAINST COMPETITORS

NOMINAL PRESSURE & RATED BURST PRESSURE OF 0.014"

| SIZE | COMPETITORS | COMPETITORS | | ENDOGLIDE | |
|----------------|---------------------|------------------|--------|------------------|--------|
| | | NOMINAL PRESSURE | RBP | NOMINAL PRESSURE | RBP |
| ø 2.0 x 220 mm | CORDIS SLEEK 2 | 5 bar | 15 bar | 9 bar | 18 bar |
| ø 1.5 x 20 mm | MEDTRONIC AMPHIRION | 7 bar | 14 bar | 9 bar | 18 bar |
| ø 3.0 x 20 mm | EV3 NANOCROSS | 10 bar | 14 bar | 9 bar | 18 bar |
| ø 4.0 x 20 mm | ABBOTT ARMADA | 8 bar | 14 bar | 9 bar | 18 bar |

ENDOGLIDE™ PRESENTS HIGHER RATED BURST PRESSURE THAN COMPETITORS PRODUCTS

ENDOGLIDE™




COMPARATIVE STUDIES
VERIFY THE SUPERIORITY OF OUR
ENDOGLIDE™ PRODUCTS

TECHNICAL SPECIFICATIONS


| | 0.014" | 0.018" | 0.035" |
|-------------------------------|---|---|---|
| MODE OF EXCHANGE | over-the-wire (OTW) | | |
| USABLE CATHETER LENGTH (UCL) | 120 cm 150 cm | 45 cm 80 cm 135 cm | 45 cm 80 cm 120 cm 150 cm |
| NOMINAL PRESSURE | 9 bar | 9 bar | 7 bar |
| RATED BURST PRESSURE | 18 bar | 18 bar | 13 or 10 bar (depending on ø) |
| BALLOON CHARACTERISTIC | semi compliant | | |
| COATING | hydrophilic coating | | |
| PLATINUM X-RAY BALLOON MARKER | 2 marker - 1x proximal and 1x distal end of balloon | | |
| RECOMMENDED INTRODUCER SHEATH | up to ø 4.5mm = 4F for ø 5mm = 5F | up to ø 4.5mm = 4F for ø 5mm to 6mm = 5F for ø 6.5mm & 7mm = 6F | up to ø 3.5mm = 5F for ø 4mm to 9mm = 6F for ø 10mm = 7F for ø 12mm = 9F |

0.014" PRODUCTS: MEAN MAXIMUM TRACKABILITY FORCE

PRODUCT VARIANT IN COMPARISON:
ø 3.0 x 20mm

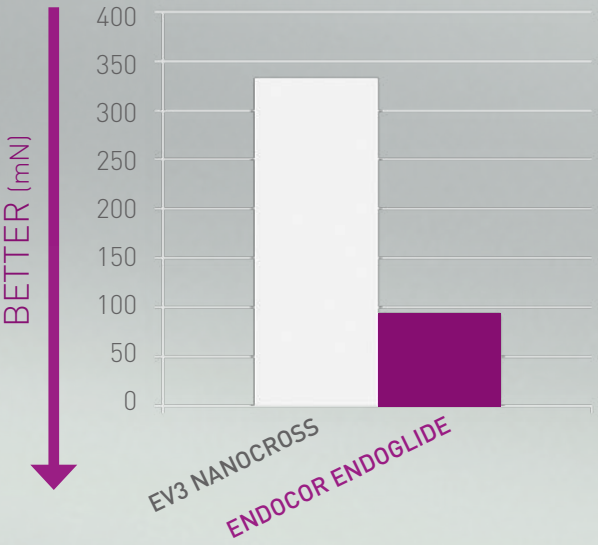


PRODUCT VARIANT IN COMPARISON:
ø 1.5 x 20mm

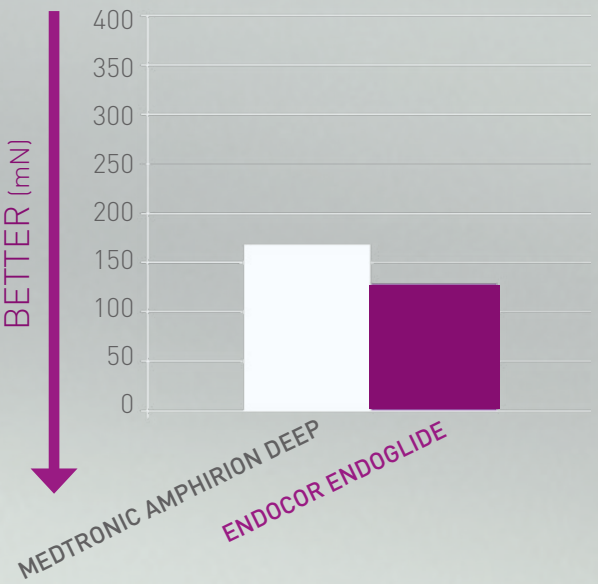

BETTER TRACKABILITY AND FASTER
DEFLATION TIME OF ENDOGLIDE™ COMPARED
TO OTHER COMPETITORS PRODUCTS


DEFLATION TIME

PRODUCT VARIANT IN COMPARISON:
0.014", ø 3.0 x 20mm

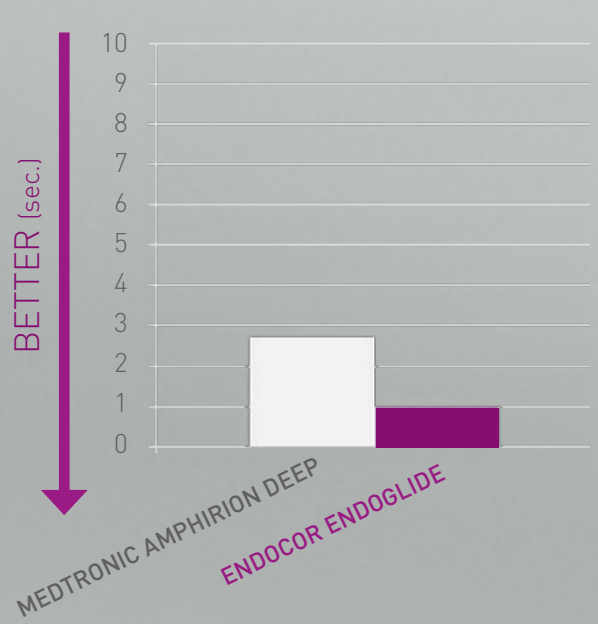

MEAN MAXIMUM TRACKABILITY FORCE

PRODUCT VARIANT IN COMPARISON:
0.014", 2.0 x 220mm

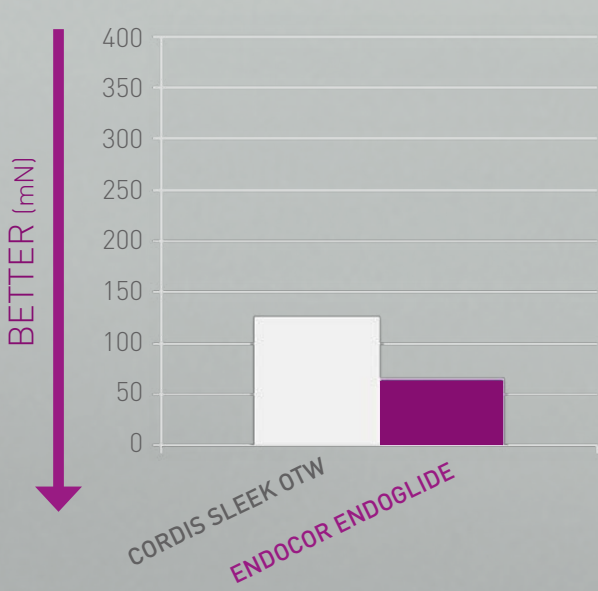

ENDOGLIDE

PRODUCT ORDERING INFORMATION: 0.035"

| PRODUCT RANGE FOR ENDOGLIDE 0.035" | | | | | | | | | | | | | | |
|---|---|----|----|----|----|----|----|-----|-----|-----|-----|-----|------|-----|
| Usable catheter length 45cm* - 80cm - 120cm - 150cm | | | | | | | | | | | | | | |
| Balloon ø (mm) | Balloon length (mm) | | | | | | | | | | | | | |
| | 10 | 15 | 20 | 30 | 40 | 60 | 80 | 100 | 120 | 150 | 180 | 200 | 250* | 280 |
| 3.00 | | | | | | | | | | | | | | |
| 3.50 | | | | | | | | | | | | | | |
| 4.00 | | | | | | | | | | | | | | |
| 4.50 | PLEASE CHECK OUR WEBSITE WWW.ENDOCOR.COM FOR DETAILED ORDERING INFORMATION | | | | | | | | | | | | | |
| 5.00 | | | | | | | | | | | | | | |
| 5.50 | | | | | | | | | | | | | | |
| 6.00 | | | | | | | | | | | | | | |
| 7.00 | | | | | | | | | | | | | | |
| 8.00 | | | | | | | | | | | | | | |
| 9.00 | | | | | | | | | | | | | | |
| 10.00 | | | | | | | | | | | | | | |
| 12.00 | | | | | | | | | | | | | | |

* minimum order quantity 25 pieces

PRODUCT ORDERING INFORMATION: 0.018"

| PRODUCT RANGE FOR ENDOGLIDE 0.018" | | | | | | | | | | | | | | |
|--|---|----|----|----|----|----|----|-----|-----|-----|-----|-----|-----|--|
| Usable catheter length 45cm - 80cm - 135cm | | | | | | | | | | | | | | |
| Balloon ø (mm) | Balloon length (mm) | | | | | | | | | | | | | |
| | 10 | 15 | 20 | 30 | 40 | 60 | 80 | 100 | 120 | 150 | 200 | 220 | 280 | |
| 2.00 | | | | | | | | | | | | | | |
| 2.50 | | | | | | | | | | | | | | |
| 3.00 | PLEASE CHECK OUR WEBSITE WWW.ENDOCOR.COM FOR DETAILED ORDERING INFORMATION | | | | | | | | | | | | | |
| 3.50 | | | | | | | | | | | | | | |
| 4.00 | | | * | | * | | | | | | | | | |
| 4.50 | | | | | | | | | | | | | | |
| 5.00 | | | * | | * | | | | | | | | | |
| 5.50 | | | | | | | | | | | | | | |
| 6.00 | | | * | | * | | | | | | | | | |
| 7.00 | | | * | | * | | | | | | | | | |

* with usable catheter length of 45 cm & 80 cm only available as ENDOGLIDE DF



Ø 1.25mm - 12mm

USABLE CATHETER LENGTHS: 45cm, 80cm, 120cm, 135cm, 150cm

VERSATILITY FOR YOUR CLINICALLY DAILY ROUTINE: ENDOGLIDE™ FEATURES AN INCREDIBLE VARIETY OF SIZES

PRODUCT ORDERING INFORMATION: 0.014"

| PRODUCT RANGE FOR ENDOGLIDE 0.014" | | | | | | | | | | | | | | | | |
|--|---|----|----|----|----|----|----|-----|-----|-----|-----|-----|-----|-----|-----|-----|
| Usable catheter length 120 cm - 150 cm | | | | | | | | | | | | | | | | |
| Balloon ø (mm) | Balloon length (mm) | | | | | | | | | | | | | | | |
| | 10 | 15 | 20 | 30 | 40 | 60 | 80 | 100 | 120 | 140 | 150 | 170 | 200 | 220 | 250 | 280 |
| 1.25 | | | * | | * | * | * | * | * | | * | * | * | * | | |
| 1.50 | | | * | | * | * | * | * | * | | * | * | * | * | | |
| 2.00 | PLEASE CHECK OUR WEBSITE WWW.ENDOCOR.COM FOR DETAILED ORDERING INFORMATION | | | | | | | | | | | | | | | |
| 2.50 | | | | | | | | | | | | | | | | |
| 3.00 | | | * | | * | * | * | * | * | | * | * | * | * | | |
| 3.50 | | | * | | * | * | * | * | * | | * | * | * | * | | |
| 4.00 | | | * | | * | * | * | * | * | | * | * | * | * | | |
| 4.50 | | | * | | * | * | * | * | * | | * | * | * | * | | |
| 5.00 | | | * | | * | * | * | * | * | | * | * | * | * | | |

* with usable catheter length of 150 cm only available as ENDOGLIDE BTK

0.035" 0.018" 0.014"

WITH AN IMPRESSIVE RANGE OF SIZES

LOWEST TIP ENTRY
PROFILE IN THE WORLD

Manufacturer:

ENDOCOR GmbH

GERMANY

info@endocor.com

www.endocor.com

Distributed by:

CE 0481

© ENDOCOR

All rights reserved. Technical specification and data are subject to modification, revision and improvement without notice. RJIS07062017