


ENDOCOR
PIONEERING INNOVATION FOR LIFE™

PHANTOM™

6F / 7F

THROMBUS REMOVING ASPIRATION CATHETER SYSTEM & KIT

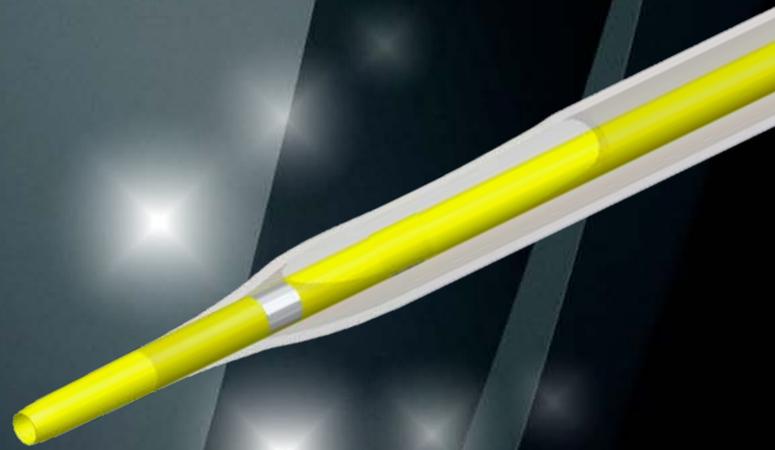

A NEW MILESTONE
IN THROMBUS REMOVAL


PHANTOM™ 6F/7F

THROMBUS REMOVING ASPIRATION CATHETER SYSTEM & KIT

This remarkable new aspiration catheter sets a new milestone in thrombus removal. Our PHANTOM™ product series complete their mission without leaving a trace.

Based on our new technology it enables the effective removal of fresh thrombus and debris in peripheral and coronary arteries by precisely directing the percutaneous suction with optimal aspiration rates to the thrombus without damaging the vessel walls.

For a wider application this innovative rapid exchange type catheter comes in the sizes 6F and 7F for treatment of both coronary and peripheral vessels.

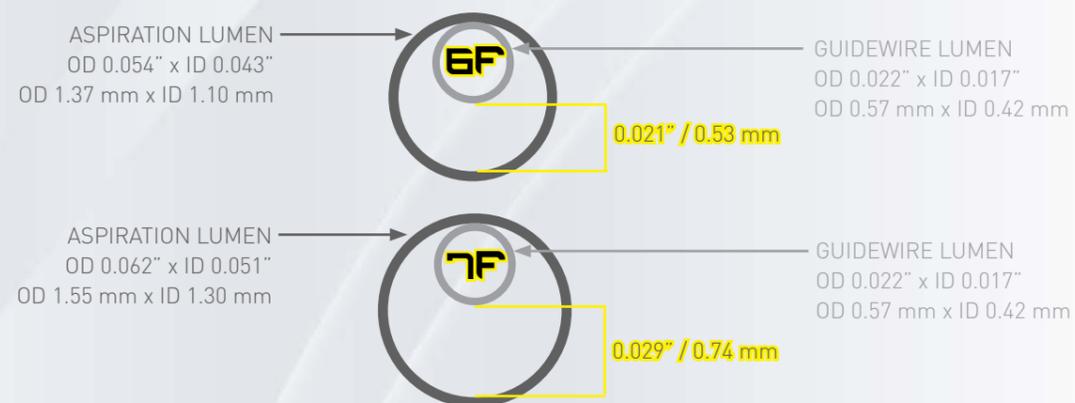


RX CATHETER SYSTEM

Supports crossing even the most tortuous paths



HIGHLIGHTS

- Improved pushability and kink resistance through pre-loaded removable stylet
- Superior flexibility and luminal volume
- Easy navigation by hydrophilic coating

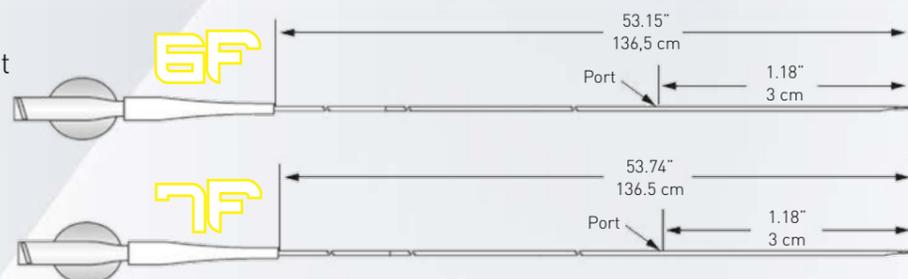


SUPERB ASPIRATION VOLUME

Our PHANTOM™ devices are designed to optimize thrombectomy: The shaft construction with a large aspiration port allows for a large extraction lumen with higher aspiration rates. This shortens treatment time and makes it more effective.

HIGHLIGHTS

- Safe procedures with excellent delivery by short tip design
- Large extraction lumen due to shaft construction, ensures optimal aspiration rates
- Effective performance by large aspiration port allows higher aspiration rates and volume

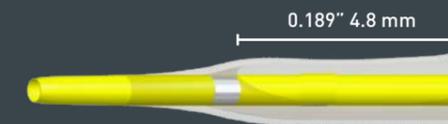

PHANTOM	6F	7F
GUIDEWIRE	0.36 mm (0.014")	0.36 mm (0.014")
RAPID EXCHANGE (RX)	30 mm (1.18")	30 mm (1.18")
WORKING LENGTH	136.5 cm (53.15")	136.5 cm (53.74")
CATHETER TYPE	NON-BRAIDED	NON-BRAIDED
MAX. OUTER DIAMETER DISTAL	1.42 mm (0.056")	1.58 mm (0.062")
INNER DIAMETER (I.D.)	1.10 mm (0.043")	1.30 mm (0.05")
ASPIRATION PORT OPENING	4.8 mm (0.189")	4.8 mm (0.189")
ENTRY TIP DIAMETER	0.40 mm < tip ϕ < 0.50 mm	0.40 mm < tip ϕ < 0.50 mm
ASPIRATION AREA (CATH. CROSS SECTION)	0.95 mm ²	1.33 mm ²
RADIOPAQUE MARKER BAND (DISTANCE TIP TO MARKER)	3.5 mm (0.14")	3.5 mm (0.14")
EXTRACTION AREA (PORT)	3,86 mm ² (length 4,8 mm, max. width 1,0 mm)	4,80 mm ² (length 4,8 mm, max. width 1,3 mm)
ORDERING INFORMATION	PM6F	PM7F

CONCEPTUALIZED FOR SAFE PROCEDURES

The short tip and atraumatic design of the PHANTOM™ series facilitate delivery. Directed suction with high aspiration rates remove debris and thrombus but do not harm the neighboring healthy tissue.

BRILLIANT TRACTABILITY,
PUSHABILITY AND ASPIRATION VOLUME.

PHANTOM 6F/7F

A NEW MILESTONE IN THROMBUS REMOVAL

Our kit includes 1x aspiration catheter, 2x aspiration syringes,
1x filter basket, 1x tubing set.

Manufacturer:

ENDOCOR GmbH

GERMANY

info@endocor.com

www.endocor.com

Distributed by:



CE 0481

© ENDOCOR

All rights reserved. Technical specification and data
are subject to modification, revision and improvement
without notice. RJIS22062019